




GC CAPILLARY COLUMN APPLICATIONS:

Flavors/Fragrances

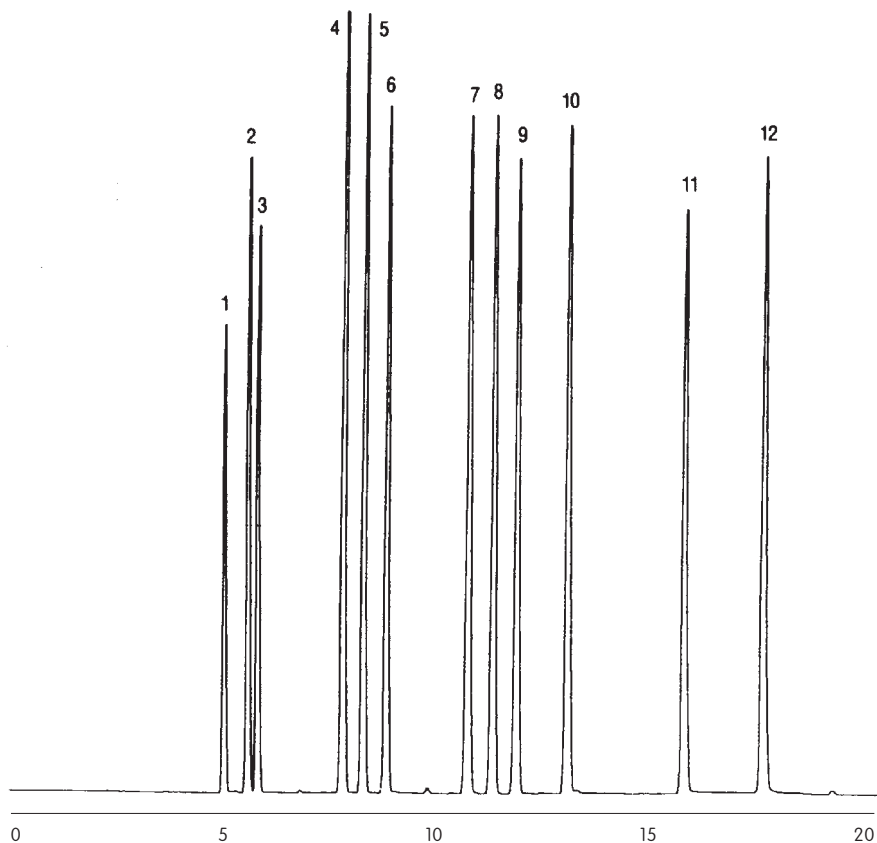

Alcohols

COLUMN: BTR-CW, Carbowax/BTR -
Oxygen resistant
30M. x 0.53mm I.D. x 1.0 μ m film
Cat. No.: BTR-CW-30V-1.0F

Temperature: 50° (5°/min.) - 100°C
Injector: 220°C
Detector: 250°C , FID
Carrier Gas: 21.7 cm/sec., Helium

- | | | |
|-----------------------|------------------------|------------------------|
| 1. methanol | 5. 2-butanol | 9. 2-pentanol |
| 2. 2-propanol | 6. 1-propanol | 10. 1-butanol |
| 3. ethanol | 7. 2-methyl-1-propanol | 11. 3-methyl-1-butanol |
| 4. 2-methyl-2-butanol | 8. 3-pentanol | 12. 1-pentanol |



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com