




GC CAPILLARY COLUMN APPLICATIONS:

Flavors/Fragrances

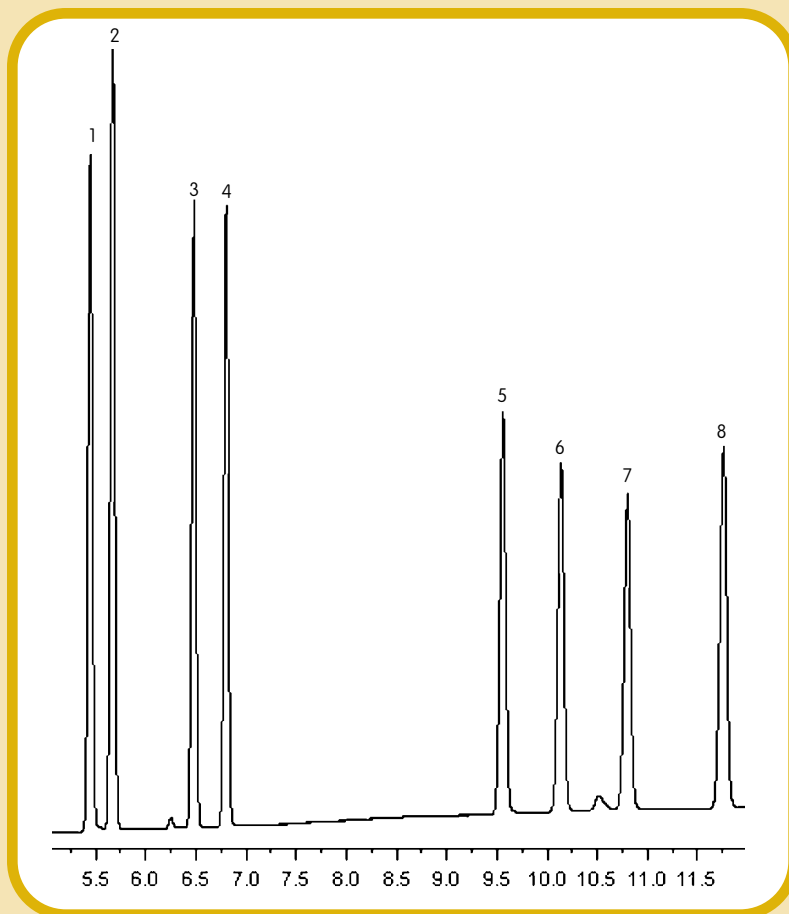

Alditol Acetates

COLUMN: 007-23, (78% Cyanopropyl) Methylpolysiloxane
30M. x 0.25mm x 0.25 μ m film
Cat. No.: 007-23-30-0.25F

Temperature: 220° (2 min. hold)
(5°/min.) - 250°C

Injector: 225°C

Detector: 250°C, FID

Carrier Gas: 18 psi, He

Injection: 0.6 μ L, split

1. Rhamnitol
2. Fucitol
3. Ribitol
4. Arabinitol
5. Mannitol
6. Galactitol
7. Glucitol
8. Inositol



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com