




GC CAPILLARY COLUMN APPLICATIONS: Steroids

Anabolic Steroids

COLUMN: 007-5*, (5% Phenyl)
Methylpolysiloxane
25M. x 0.53mm I.D. x 1.0µm film
Cat. No.: 007-5-25V-1.0F

Temperature: 260°C
Injector: 250°C
Detector: 310°C, FID
Carrier Gas: 7 ml/min., Helium

- | | |
|-------------|-----------------|
| 1. dianabol | 4. testosterone |
| 2. anavar | 5. dolderone |
| 3. diandron | |

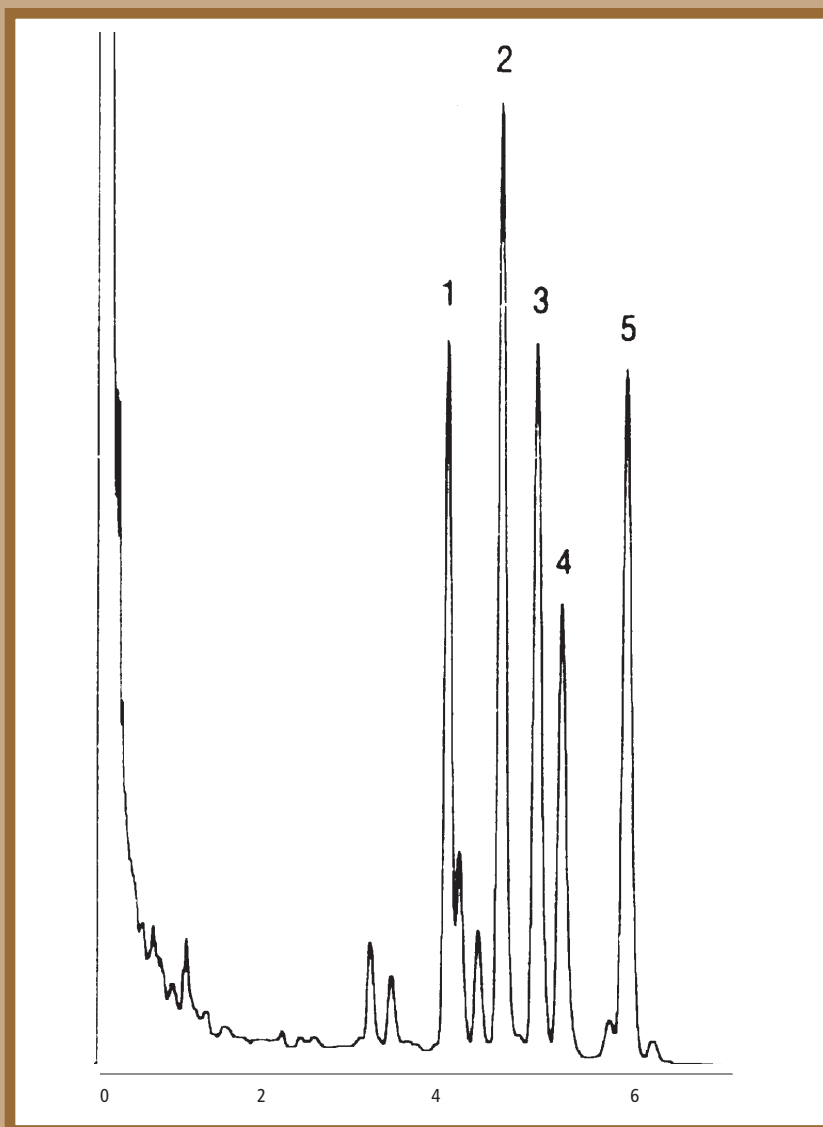

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com