




# GC CAPILLARY COLUMN APPLICATIONS: Environmental

**Aroclor  
1254**

**COLUMN:** 007-608, Cyanopropyl Phenyl  
Methylpolysiloxane  
30M. x 0.53mm I.D. x 0.8 $\mu$ m film  
Cat. No.: 007-608-30V-0.8F

**Temperature:** 140° (8°/min.) - 260°C  
**Injector:** 240°C  
**Detector:** 300°C, ECD  
**Carrier Gas:** 10.5 ml/min., Helium

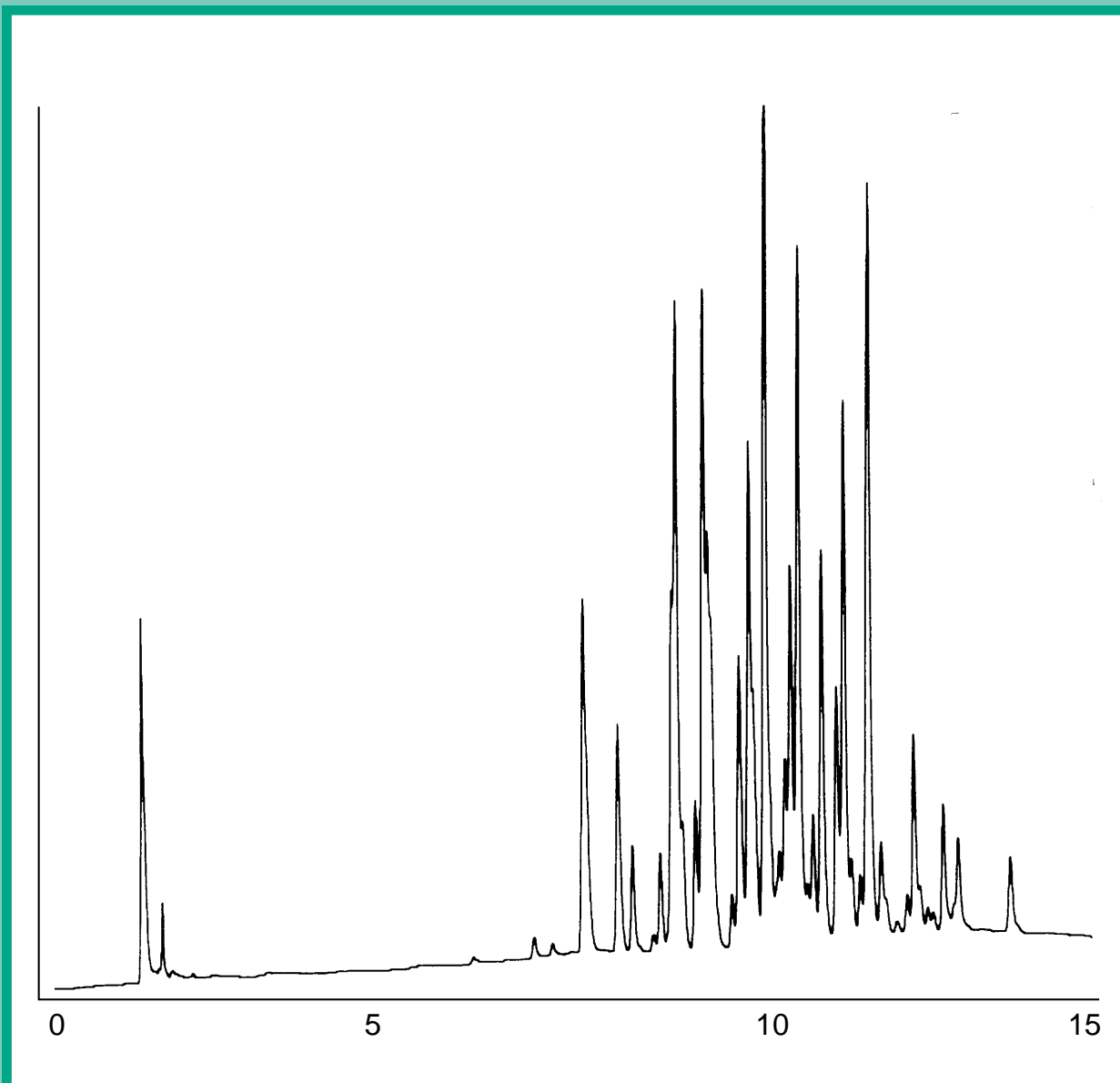

**QUADREX**



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com) • WEB: [www.quadrexcorp.com](http://www.quadrexcorp.com)