




GC CAPILLARY COLUMN APPLICATIONS: Triglycerides

Butter Triglycerides

COLUMN: 400-65HT, (65% Phenyl)
Methylpolysiloxane - Aluminum-Clad
15M. x 0.25mm I.D. x 0.1µm film
Cat. No.: 400-65HT-15-0.1F

Temperature: 250° (5°/min.) - 365°C
Injector: 370°C
Detector: 400°C, FID
Carrier Gas: 35 cm/sec., Helium

- | | | | |
|----------------|--------|--------|---------|
| 1. cholesterol | 4. ppo | 7. poo | 10. ooo |
| 2. ppp | 5. pss | 8. ssp | |
| 3. pps | 6. pso | 9. soo | |

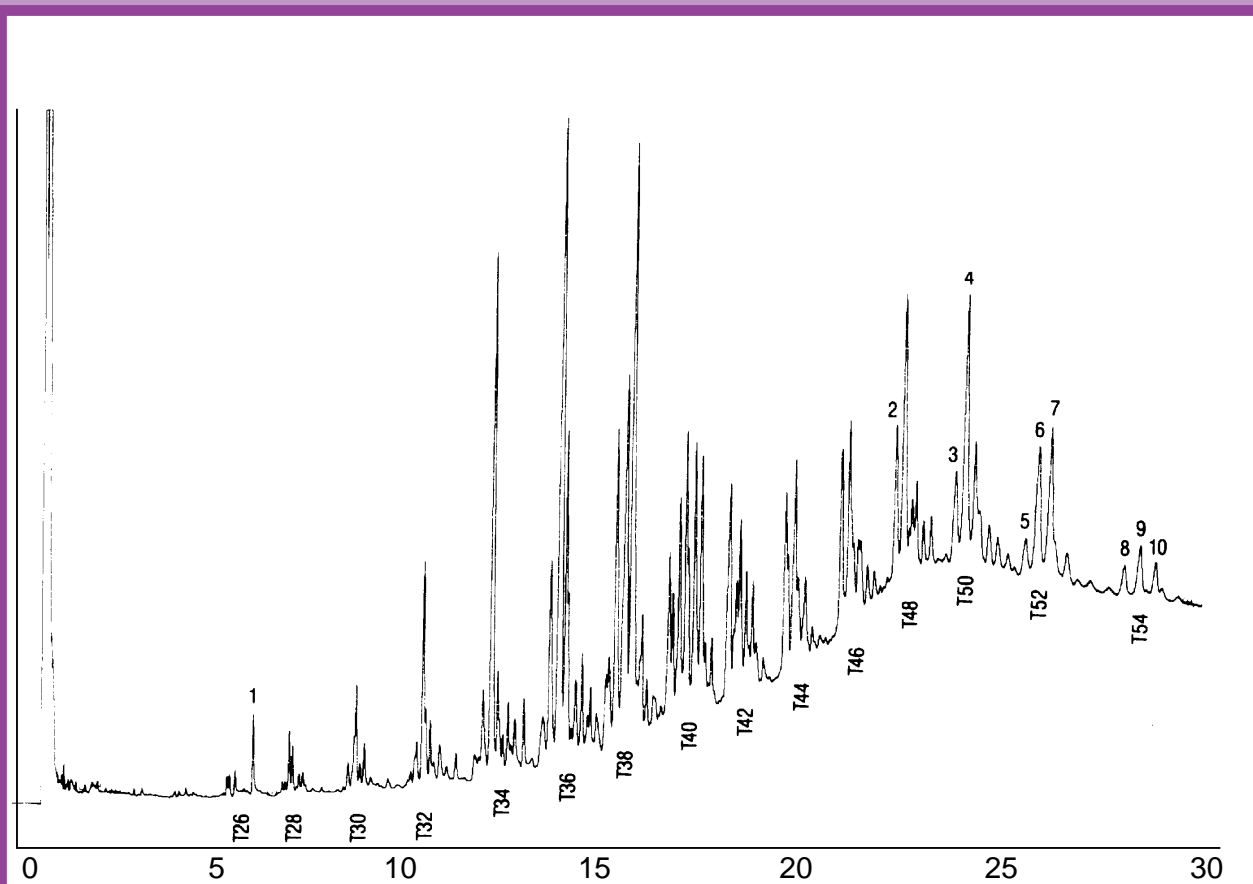

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com