




# GC CAPILLARY COLUMN

## APPLICATIONS: Ultra Alloy™ Stainless Steel

*C12 to C60  
Hydrocarbons*

**COLUMN:** UAC-DX30, Carborane siloxane  
15M. x 0.25mm I.D. x 0.15µm film  
Cat. No.: UAC-DX30-15-0.15F

**Temperature:** 100° (20°/min.) - 450°  
**Injector:** 350°C  
**Detector:** 420°C, FID  
**Carrier Gas:** 30 cm/sec., Helium

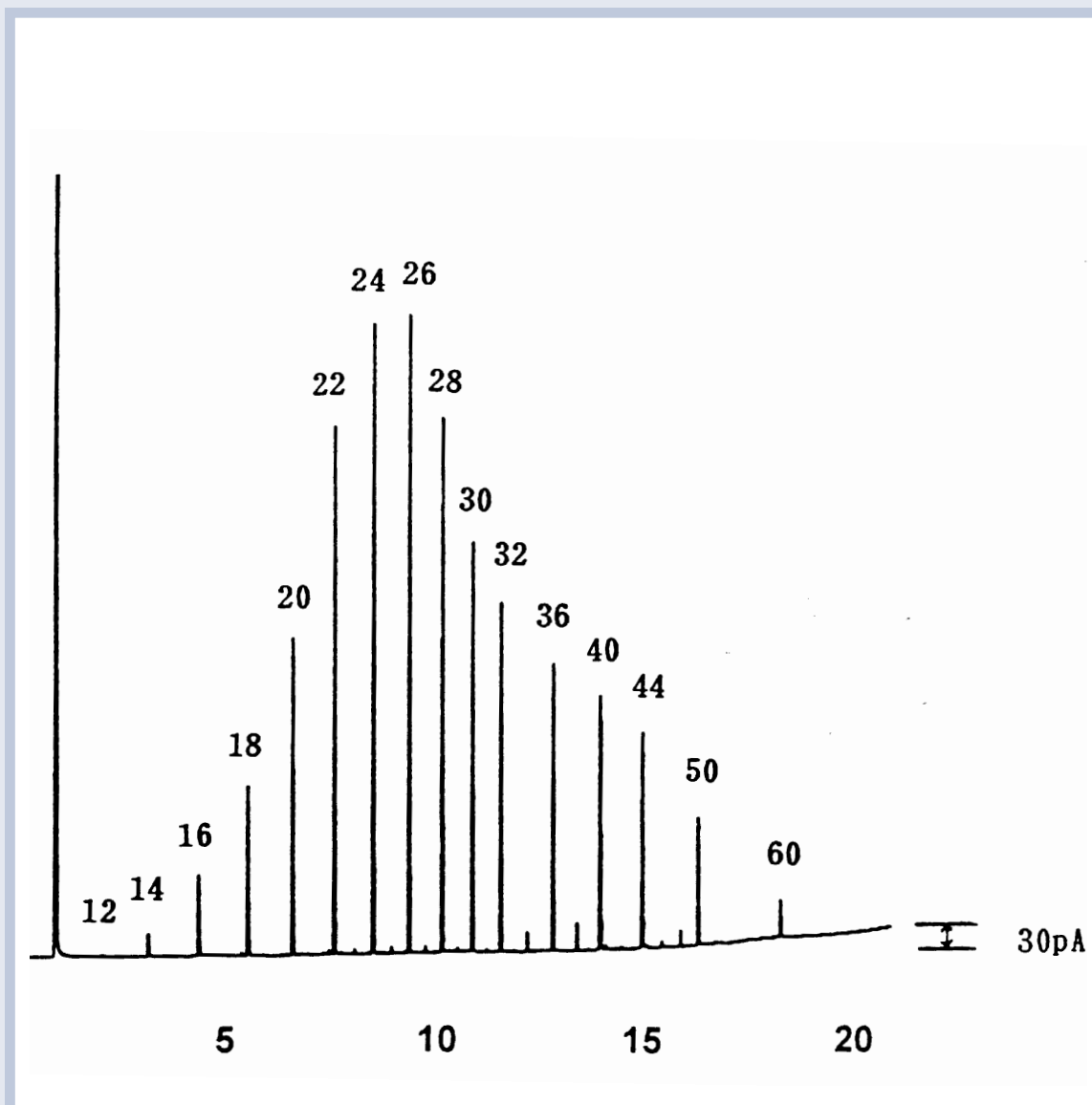

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com

**QUADREX**

