




# GC CAPILLARY COLUMN APPLICATIONS:

## Flavors/Fragrances

### Celery Seed Oil

**COLUMN:** BTR-CW, Carbowax/BTR -  
Oxygen Resistant  
30M. x 0.25mm I.D. x 0.25 $\mu$ m film  
Cat. No.: BTR-CW-30-0.25F

**Temperature:** 60° (5°/min.) - 240°C  
**Injector:** 225°C  
**Detector:** 250°C, FID  
**Carrier Gas:** 27 cm/sec., Helium

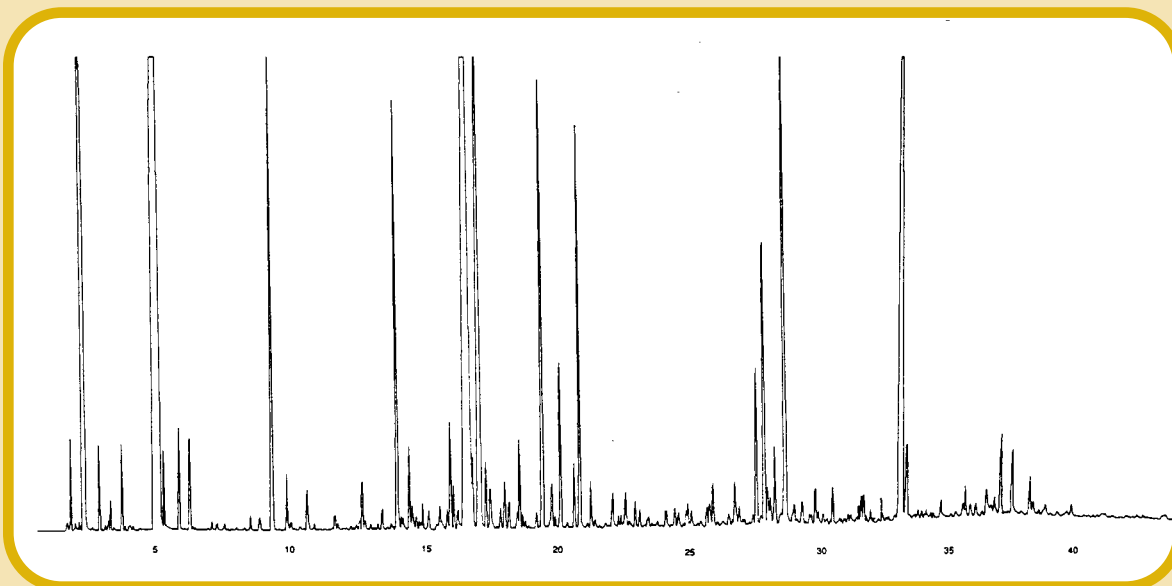

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com) • WEB: [www.quadrexcorp.com](http://www.quadrexcorp.com)