




# GC CAPILLARY COLUMN APPLICATIONS: Triglycerides

## Cocoa Butter

**COLUMN:** 400-65HT, Bonded (65% Phenyl) Methylpolysiloxane - Aluminum-Clad 15M. x 0.25mm I.D. x 0.1µm film  
Cat. No.: 400-65HT-15-0.1F

**Temperature:** 350° (0.5min. hold) (1°/min.) - 365°C  
**Injector:** 400°C  
**Detector:** 400°C, FID  
**Carrier Gas:** 35cm/sec., Helium

- |        |        |         |         |
|--------|--------|---------|---------|
| 1. ppp | 5. plp | 9. pls  | 13. soa |
| 2. mop | 6. pss | 10. sss |         |
| 3. pps | 7. pos | 11. sos |         |
| 4. pop | 8. poo | 12. soo |         |

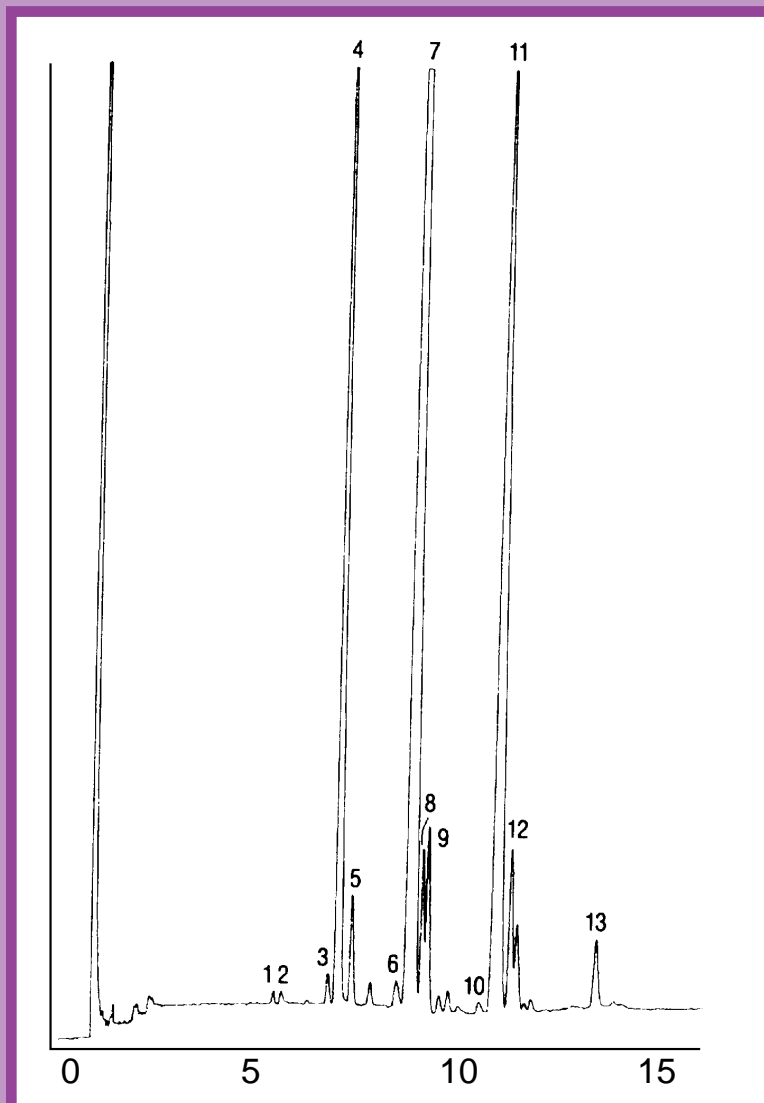

**QUADREX**



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com