




GC CAPILLARY COLUMN APPLICATIONS: Steroids

Corticosteroids

COLUMN: 007-5* (5% Phenyl)
Methylpolysiloxane
15M. x 0.53mm I.D. x 1.0 μ m film
Cat. No.: 007-5-15V-1.0F

Temperature: 150° (6°/min.) - 280°C
Injector: 250°C
Detector: 300°C, FID
Carrier Gas: 5 ml/min., Helium

- | | |
|-------------------|---------------------------------|
| 1. cortisone | 4. fludrocortisone |
| 2. prednisone | 5. 6 α -methylprednisone |
| 3. hydrocortisone | |

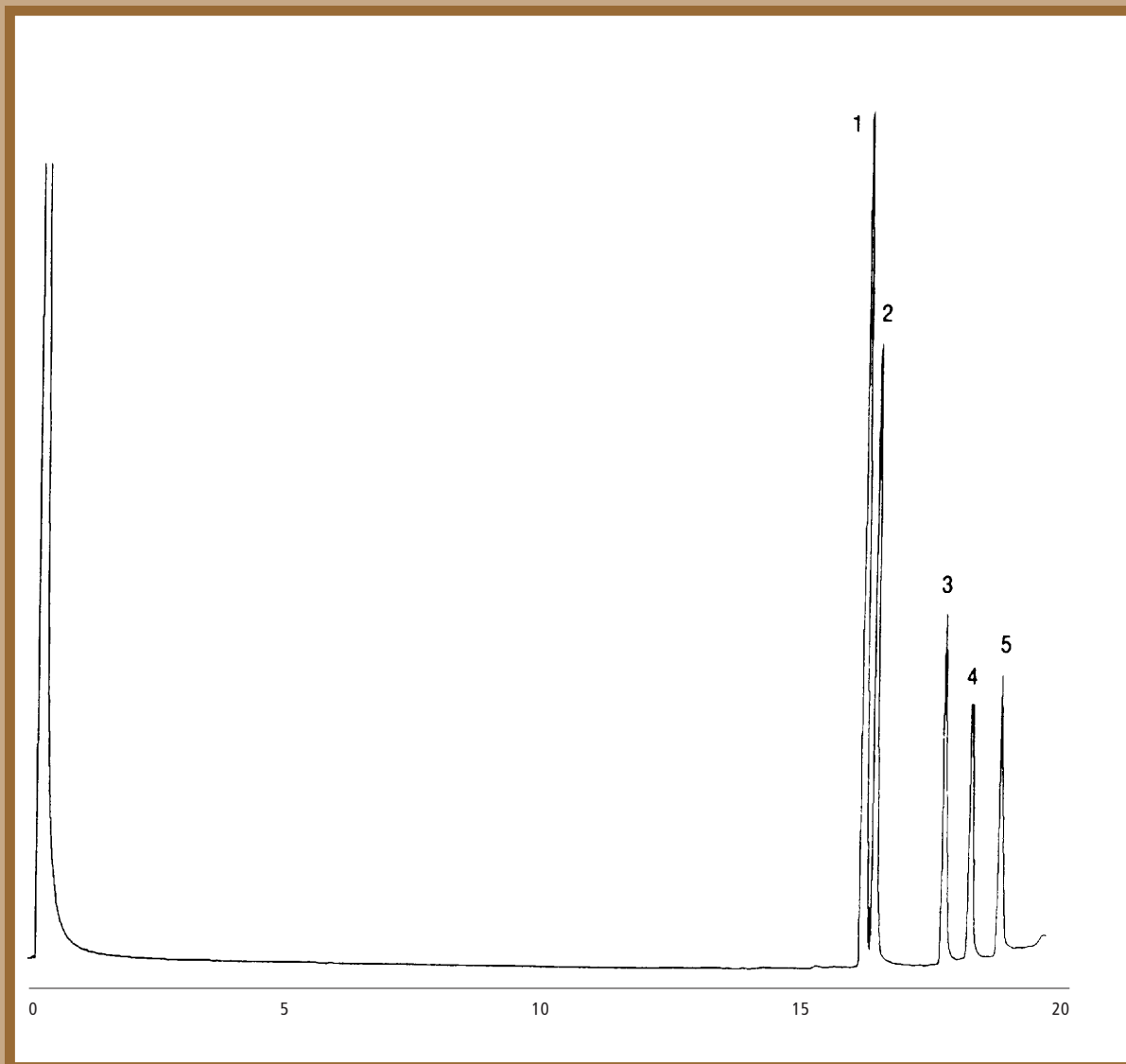

QUADREX

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com

