




GC CAPILLARY COLUMN APPLICATIONS: Environmental

DDT Isomers

COLUMN: 007-608, Cyanopropyl Phenyl Methylpolysiloxane
30M. x 0.53mm I.D. x 0.8µm film
Cat. No.: 007-608-30V-0.8F

Temperature: 140° (6°/min.) -260°C
Injector: 240°C
Detector: 300°C , ECD
Carrier Gas: 5 ml/min., Helium

- | | | |
|-------------|-------------|-------------|
| 1. o,p'-DDE | 4. m,p'-DDD | 7. p,p'-DDT |
| 2. p,p'-DDE | 5. p,p'-DDD | |
| 3. o,p'-DDD | 6. o,p'-DDT | |

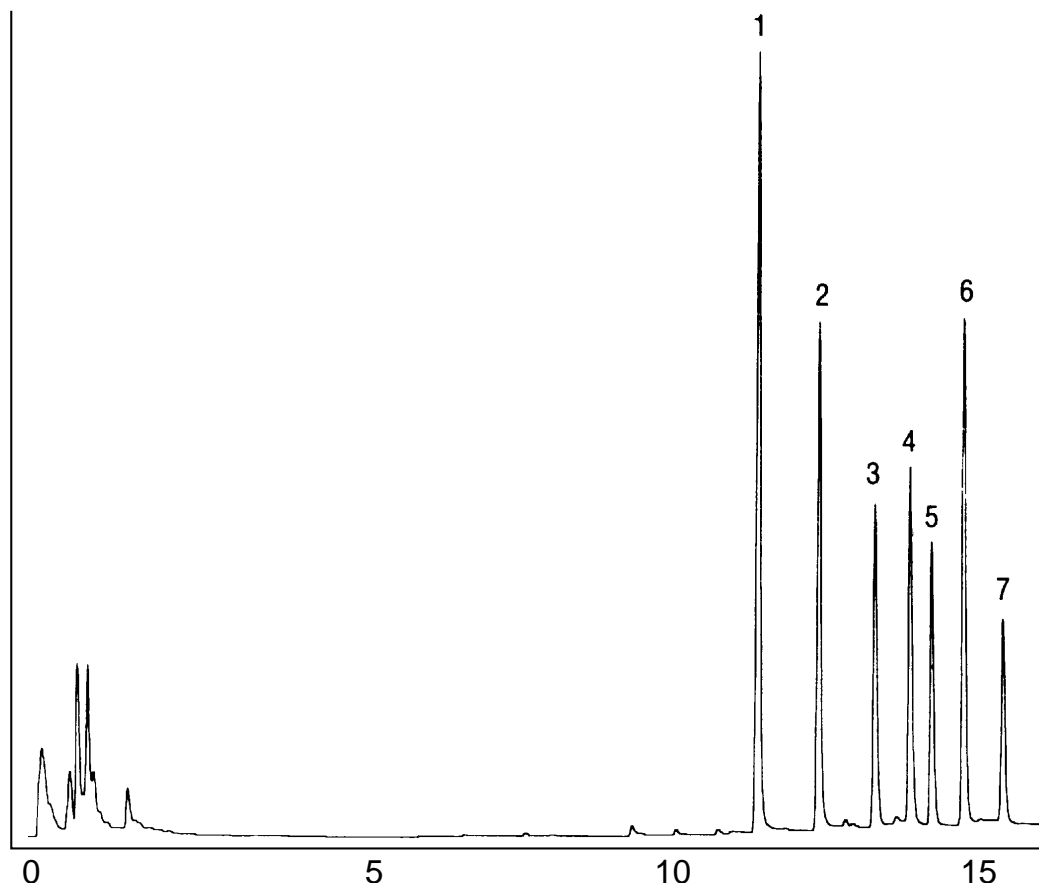

QUADREX

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com

