



GC CAPILLARY COLUMN APPLICATIONS: Environmental

Dioxins/ Furans

COLUMN:

007-23, Cyanopropyl Polysiloxane
50M. x 0.25mm I.D. x 0.25µm film
Cat. No.: 007-23-50-0.25F

Temperature:

70° (20°/min.) - 170° (3°/min.) -270°C

Injector:

230°C

Detector:


MSD

Carrier Gas:

18.5 psi., Helium

1. 1,2,3,4-TCDF
2. 2,3,7,8-TCDD

3. 1,2,3,4-TCDD
4. 2,3,7,8-TCDF

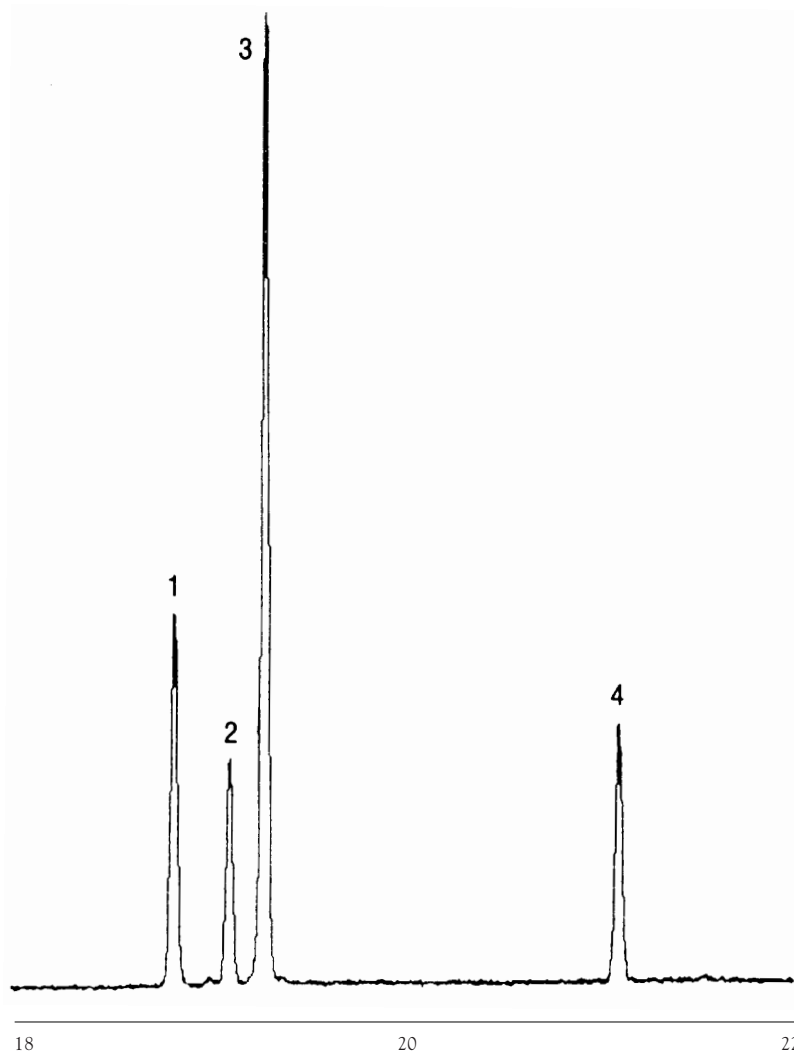

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com