




GC CAPILLARY COLUMN APPLICATIONS: Environmental

EPA Method 525 Pesticides

Column: 007-5,* (5% Phenyl) Methylpolysiloxane
30M. x 0.25mm I.D. x 0.25µm film
Cat. No.: 007-5-30-0.25F

Temperature: 60° (6°/min.) - 260°C
Injector: 230°C
Detector: MSD
Carrier Gas: 28 cm/sec., Helium

- | | |
|---------------|-----------------------|
| 1. simazine | 7. heptachlor epoxide |
| 2. atrazine | 8. g-chlordane |
| 3. lindane | 9. a-chlordane |
| 4. alachlor | 10. trans-nonachlor |
| 5. heptachlor | 11. endrin |
| 6. aldrin | 12. methoxychlor |

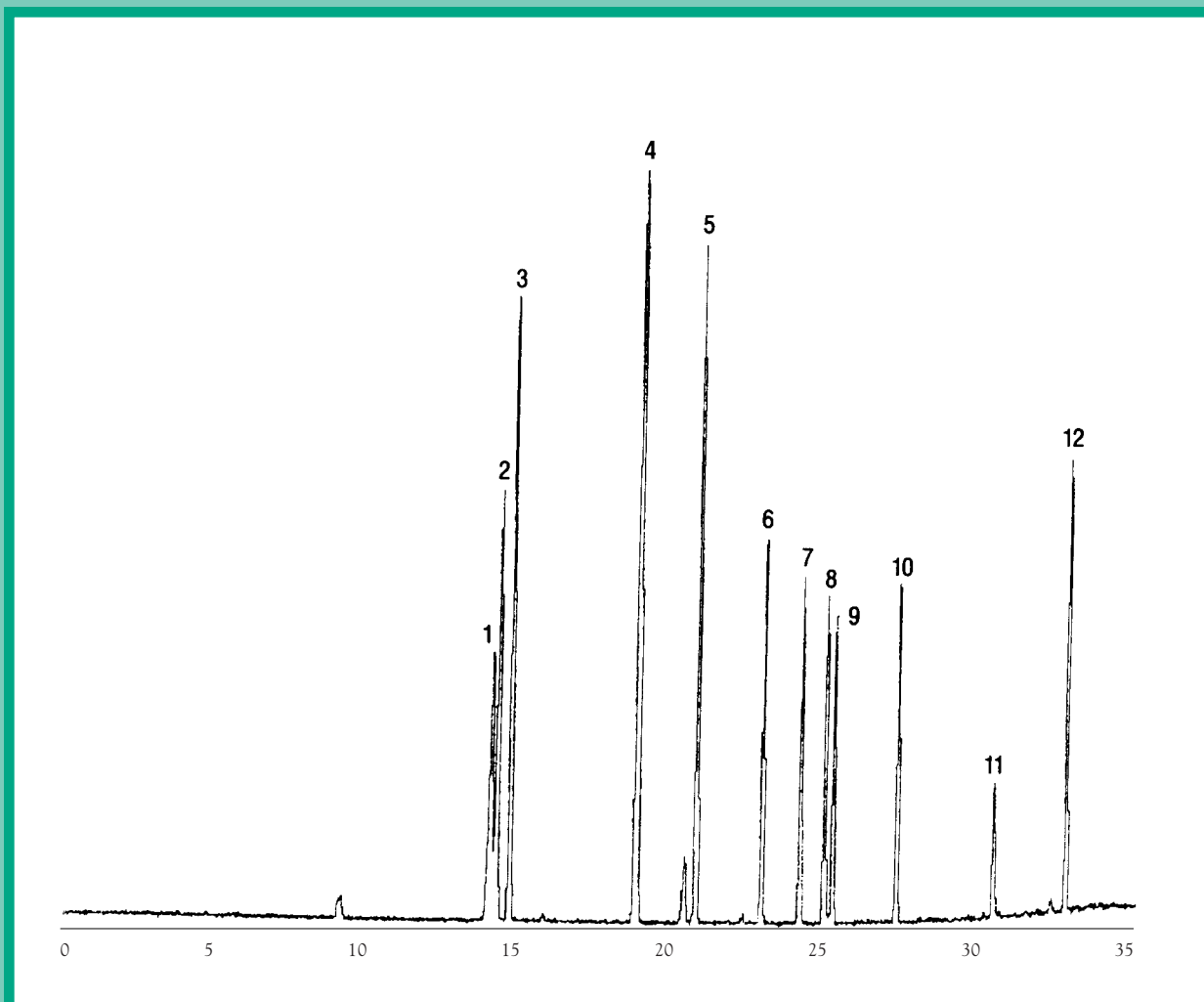

QUADREX

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com

