



GC CAPILLARY COLUMN APPLICATIONS: Environmental


EPA Method 603

COLUMN: 007-CW, Bonded Carbowax
30M. x 0.53mm I.D. x 1.0µm film
Cat. No.: 007-CW-30V-1.0F

Temperature: 50°C Isothermal
Injector: 210°C
Detector: 250°C, FID
Carrier Gas: 10 ml/min., Helium

1. acrolein

2. acrylonitrile

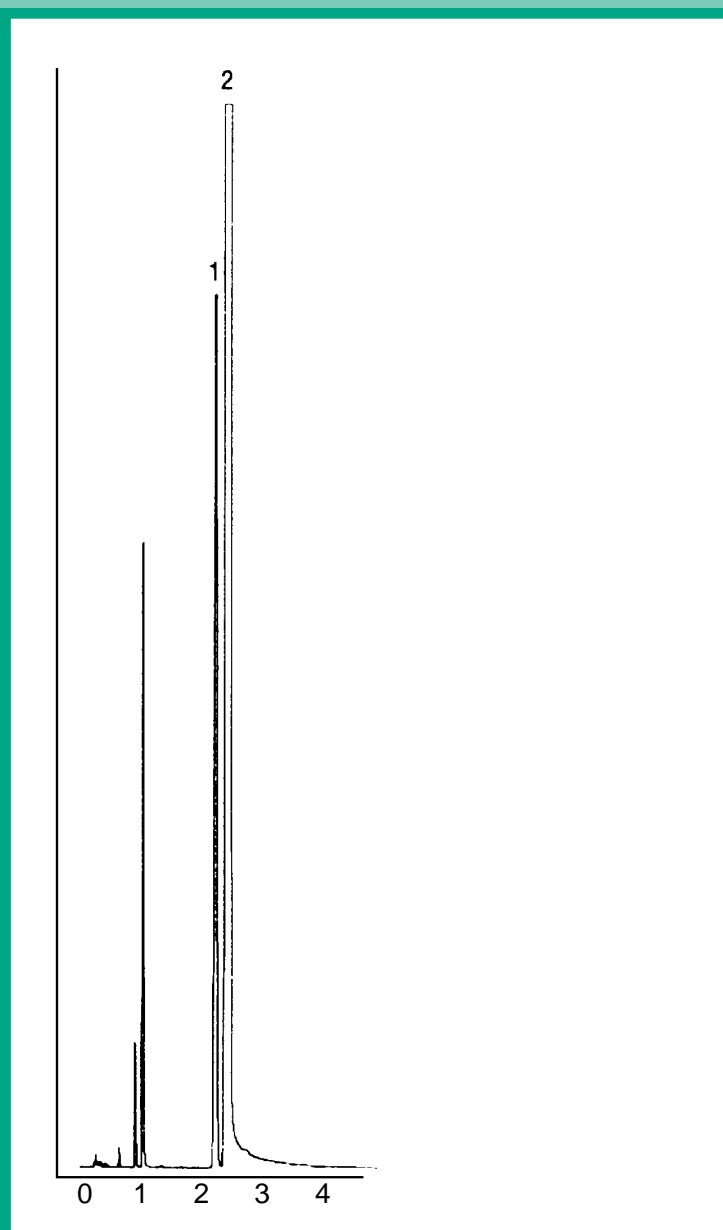

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com