




GC CAPILLARY COLUMN APPLICATIONS:

Flavors/Fragrances

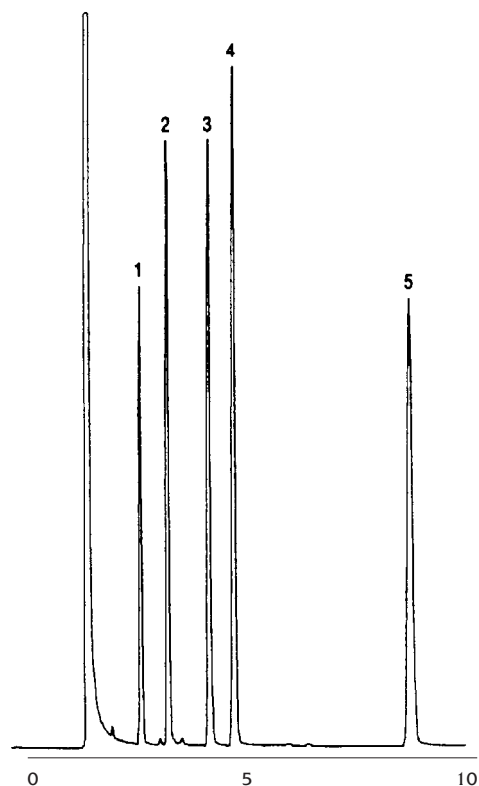

Free Fatty Acids (in Methanol)

COLUMN: 007-FFAP, Bonded FFAP
15M. x 0.32mm I.D. x 0.25 μ m film
Cat. No.: 007-FFAP-15W-0.25F

Temperature: 140°C
Injector: 220°C
Detector: 250°C, FID
Carrier Gas: 25 cm/sec., Helium

1. acetic acid
2. propionic acid
3. butyric acid
4. isovaleric acid
5. hexanoic acid