




GC CAPILLARY COLUMN APPLICATIONS: Fames

Fame Mix

COLUMN: 007-23, Bonded Cyanopropyl Polysiloxane
25M. x 0.25mm I.D. x 0.25µm film
Cat. No.: 007-23-25-0.25F

Temperature: 100° (8°/min.) - 240°C
Injector: 230°C
Detector: 300°C, FID
Carrier Gas: 28 cm/sec., Helium

1. C-14:0	5. C-16:1 cis	9. C-18:2	13. C-20:2
2. C-14:1	6. C-18:0	10. C-20:0	14. C-20:3
3. C-16:0	7. C-18:1 trans	11. C-18:3	
4. C-16:1 trans	8. C-18:1 cis	12. C-20:1	

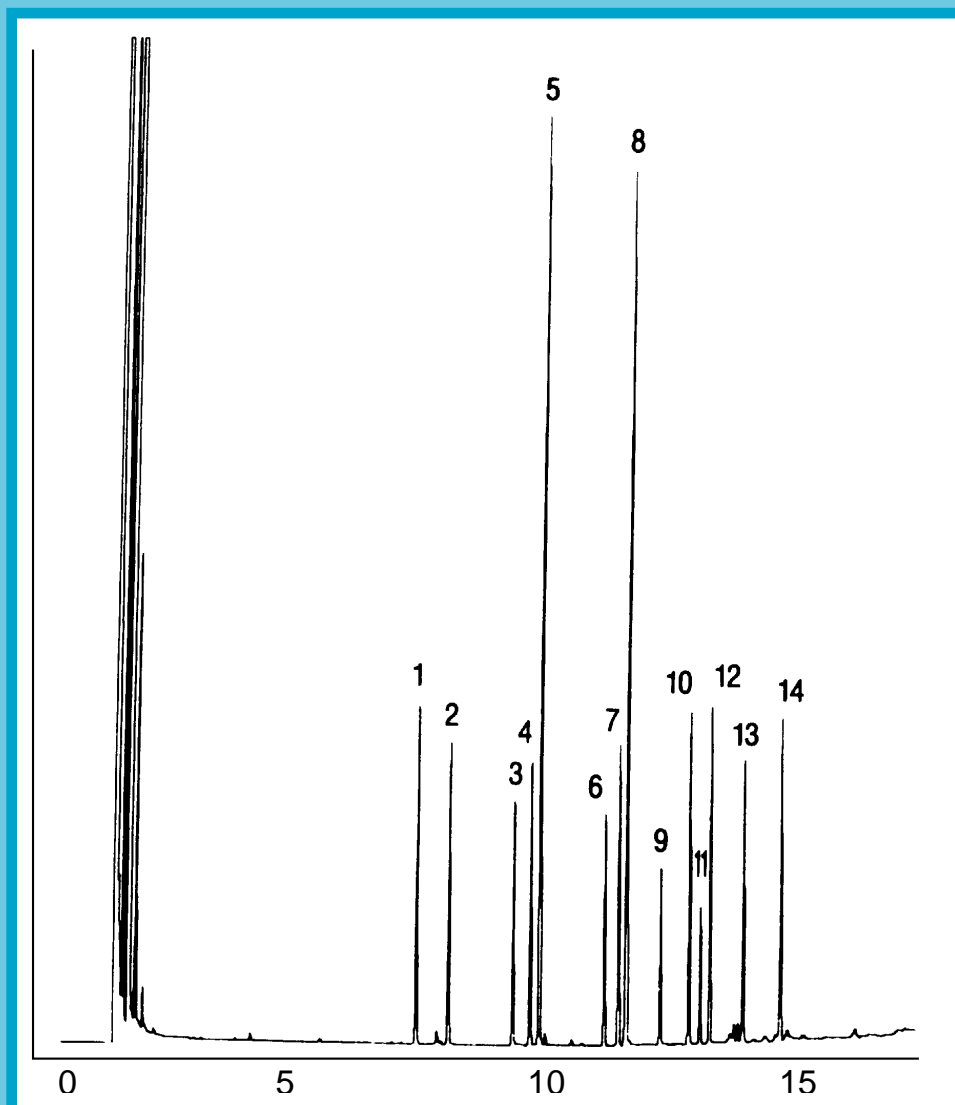

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com