



# GC CAPILLARY COLUMN

## APPLICATIONS: Ultra Alloy™ Stainless Steel

### Grob Test Mix

#### COLUMN:

UAC-CW, Carbowax  
15M. x 0.25mm I.D. x 0.25µm film  
Cat. No.: UAC-CW-15-0.25F

#### Temperature:

60° (3 min. hold) (5°/min.) - 200°C

#### Injector:

150°C (on-column)

#### Detector:

300°C, FID


#### Carrier Gas:

20 ml/min., Helium, U=66cm/sec.

1. dodecane
2. butanediol
3. octanol
4. methyl decanoate

5. dicyclohexylamine
6. methylundecanoate
7. methyl dodecanoate
8. 2,6-dimethylaniline

9. 2,6-dimethylphenol
10. ethylhexanoic acid

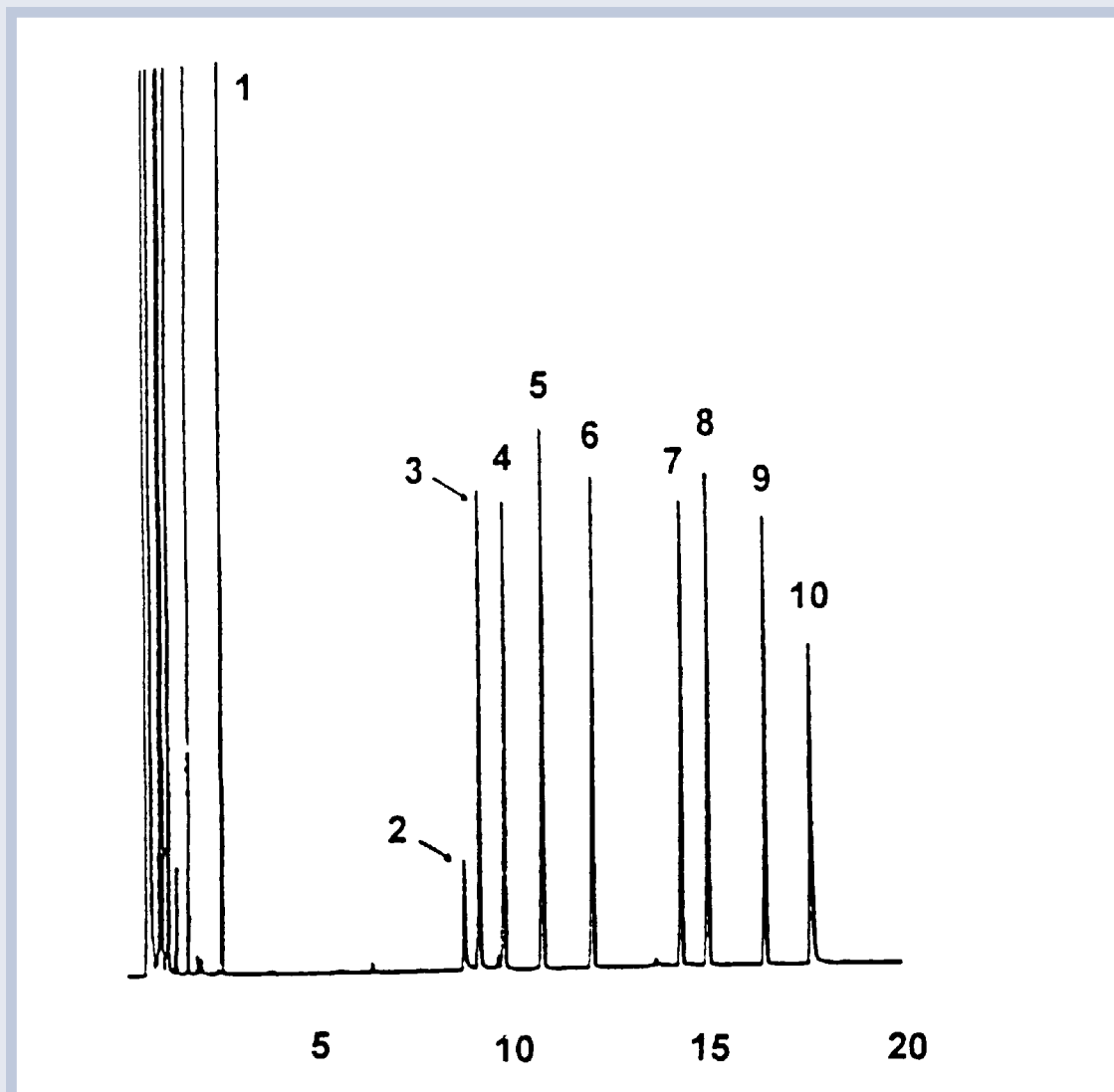

**QUADREX**



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com