



GC CAPILLARY COLUMN


APPLICATIONS: High-Temp Aluminum Clad

*Iranian
Crude
Oil*

COLUMN: 400-5HT* (5% Phenyl)
Methylpolysiloxane - Aluminum-Clad
10M. x 0.25mm I.D. x 0.1 μ m film
Cat. No.: 400-5HT-10-0.1F

Temperature: 60° (6 /min.) - 400°C
Injector: 60° - 400°C, PTV Injector, 2 μ L
Detector: 400°C, FID
Carrier Gas: 22 cm/sec., Helium

- | | |
|-------------|---------|
| 1. pristane | 4. C-20 |
| 2. phytane | 5. C-40 |
| 3. C-20 | |

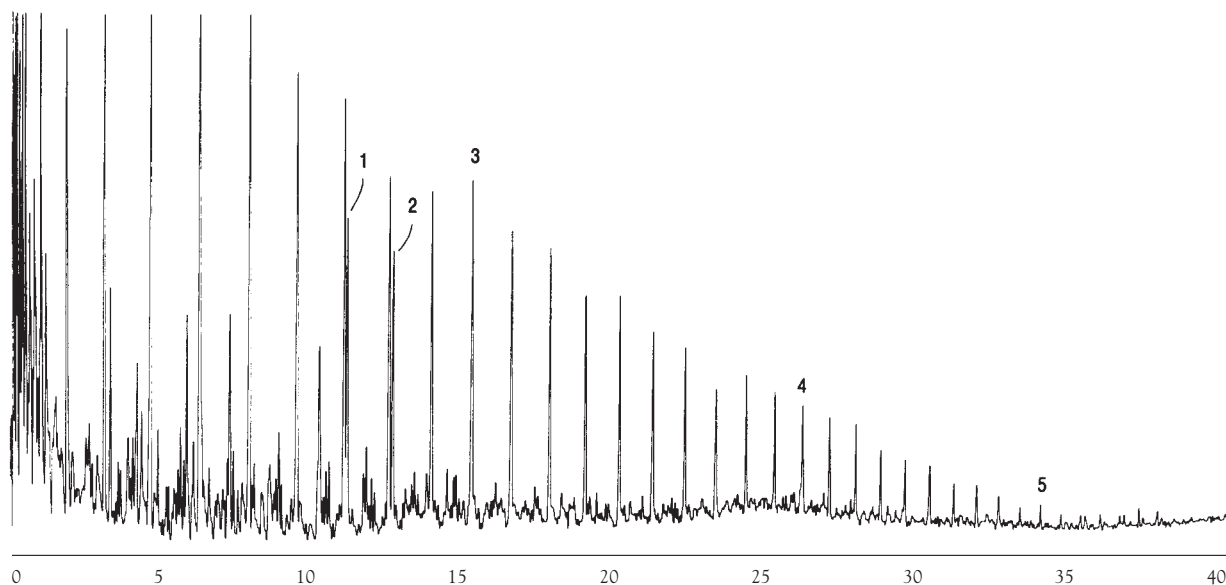

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com