




# GC CAPILLARY COLUMN APPLICATIONS: Industrial Solvents

*Isoparaffin  
Solvents  
Isopar E*

**COLUMN:** 007-1, Dimethylpolysiloxane  
30M. x 0.25mm x 0.25 $\mu$ m film  
Cat. No.: 007-1-30-0.25F

Temperature: 40° (2 min. hold) (10°/min.) - 300°C  
Injector: 225°C  
Detector: MSD, Interface temp.: 300°C  
Flow: 28cm/sec., He

ISOPAR is a trademark of Exxon Corporation.

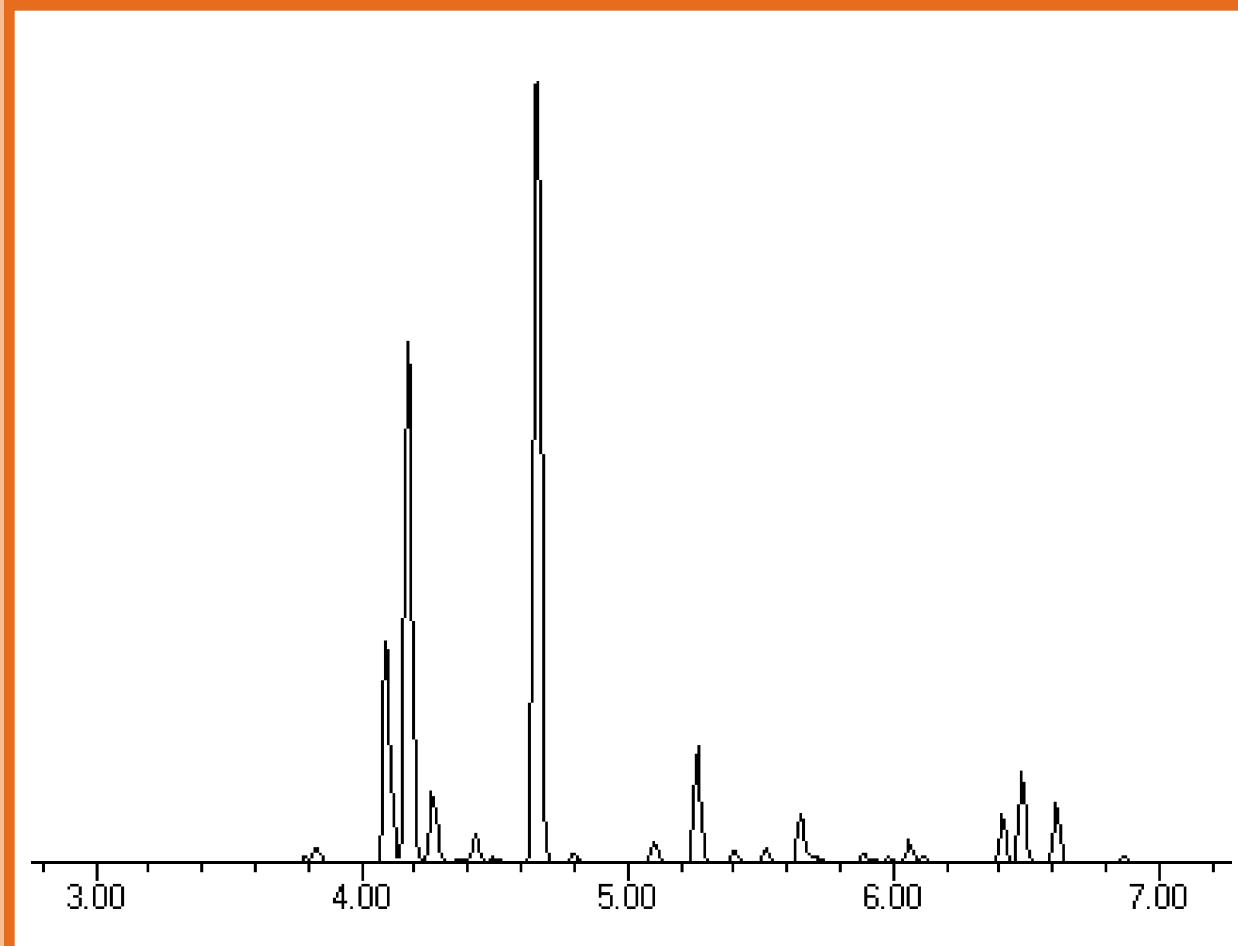

**QUADREX**



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com) • WEB: [www.quadrexcorp.com](http://www.quadrexcorp.com)