



GC CAPILLARY COLUMN

APPLICATIONS: Industrial Solvents


*Isoparaffin
Solvents*

Isopar L

COLUMN: 007-1, Dimethylpolysiloxane
30M. x 0.25mm x 0.25 μ m film
Cat. No.: 007-1-30-0.25F

Temperature: 60° (8°/min.) - 300°C
Injector: 225°C
Detector: MSD, Interface temp.: 300°C
Flow: 28cm/sec., He

ISOPAR is a trademark of Exxon Corporation.

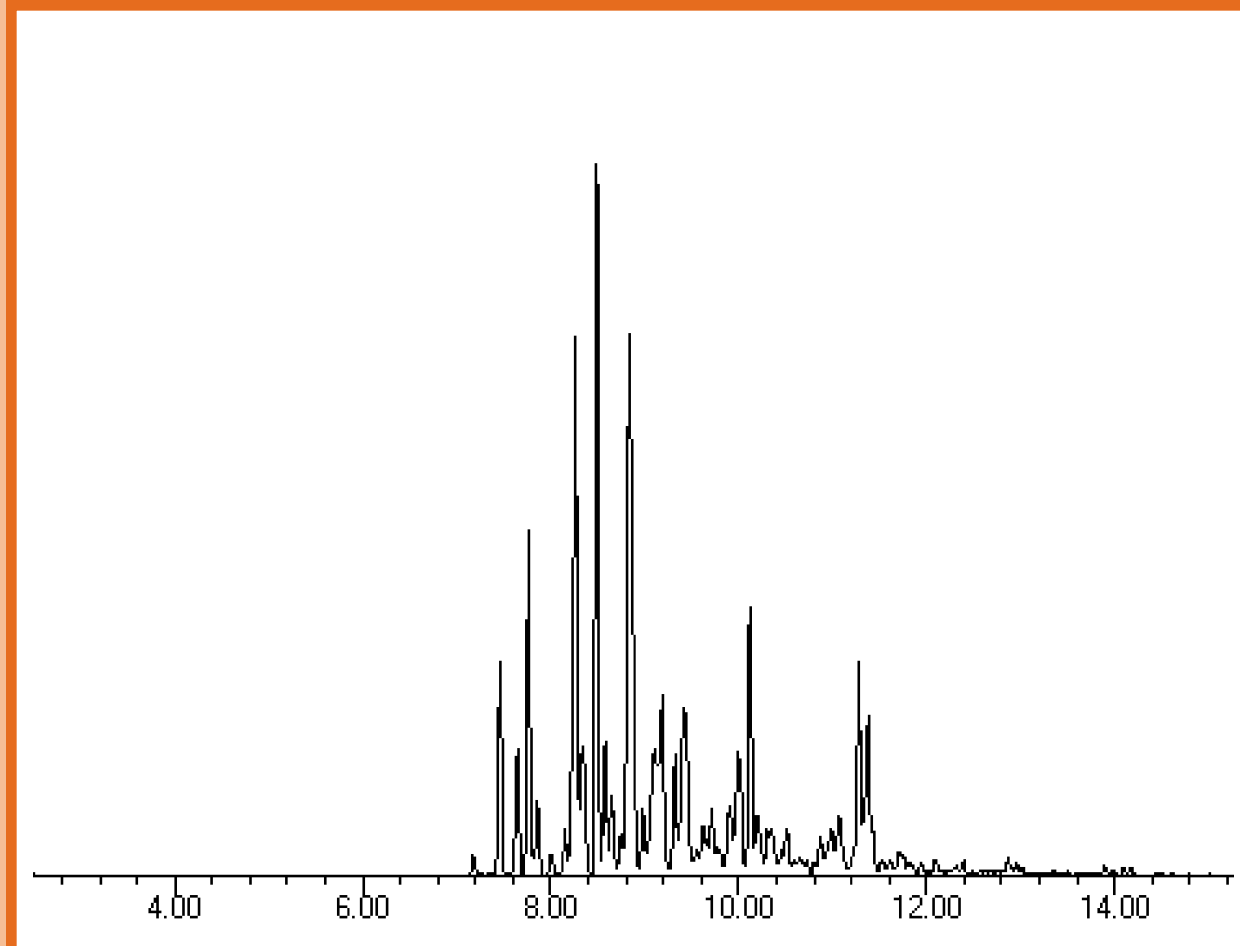

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com