




GC CAPILLARY COLUMN APPLICATIONS: Fames

Methyl Esters of Sunflower Oil

COLUMN: 007-CW, Bonded Carbowax
25M. x 0.25mm I.D. x 0.25 μ m film
Cat. No.: 007-CW-25-0.25F

Temperature: 160° (20°/min.) - 200°C
Injector: 230°C
Detector: 280°C, FID
Carrier Gas: 30 cm/sec., Helium

1. C-14:0 2. C-16:0 3. C-18:0 4. C-18:1 5. C-18:2

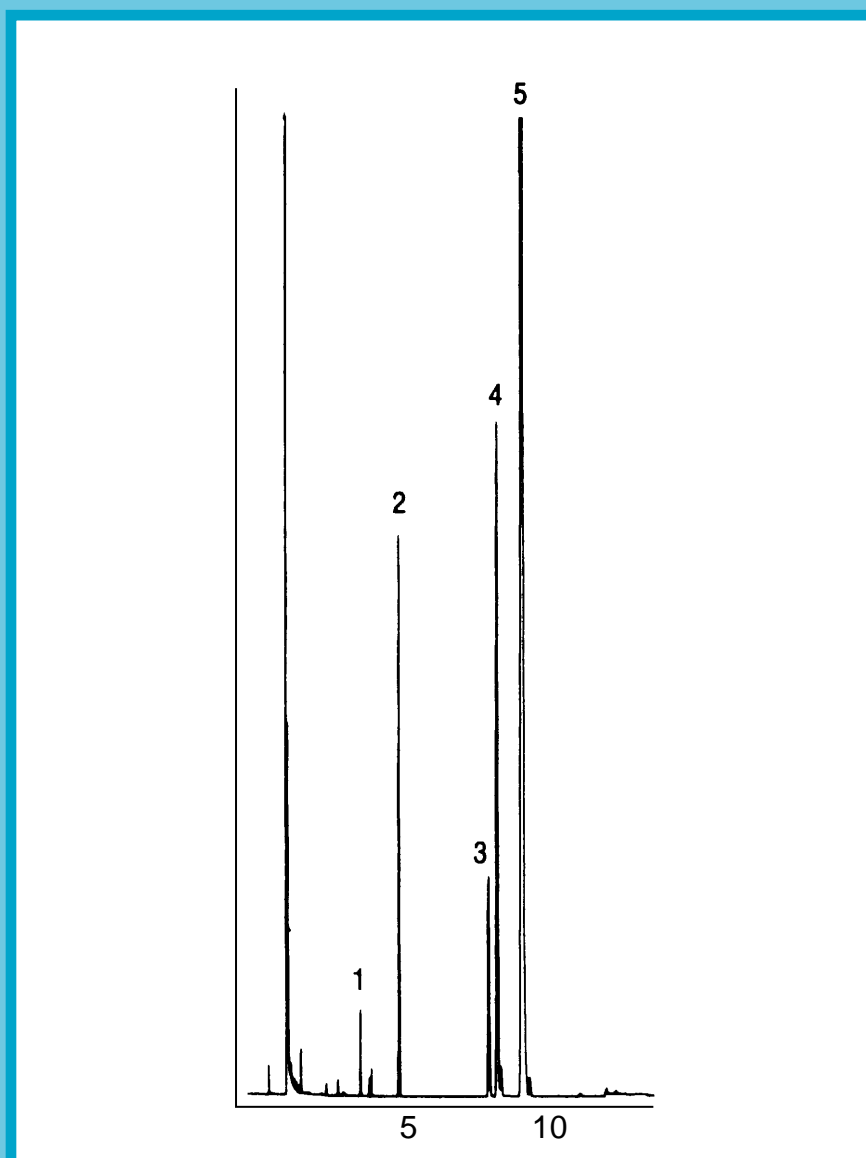

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com