



GC CAPILLARY COLUMN APPLICATIONS: **Fames**

Monoglycerides

*TMS
Derivatives*

COLUMN: 007-23, Bonded Cyanopropyl Siloxane
8.5M. x 0.25mm I.D. x 0.18 μ m film


Temperature: 80° (20°/min.) - 170°C -
(8°/min.) - 230°C

Injector: 225°C

Detector: FID

Carrier Gas: 34cm/sec., Helium

- | | |
|--------------------------|-------------------------------|
| 1. monomyristin glycerol | 5. monoleonolenin glycerol |
| 2. monopalmitin glycerol | 6. monoeicosanin glycerol |
| 3. monostearin glycerol | 7. monoeicosatrienin glycerol |
| 4. monoolein glycerol | |

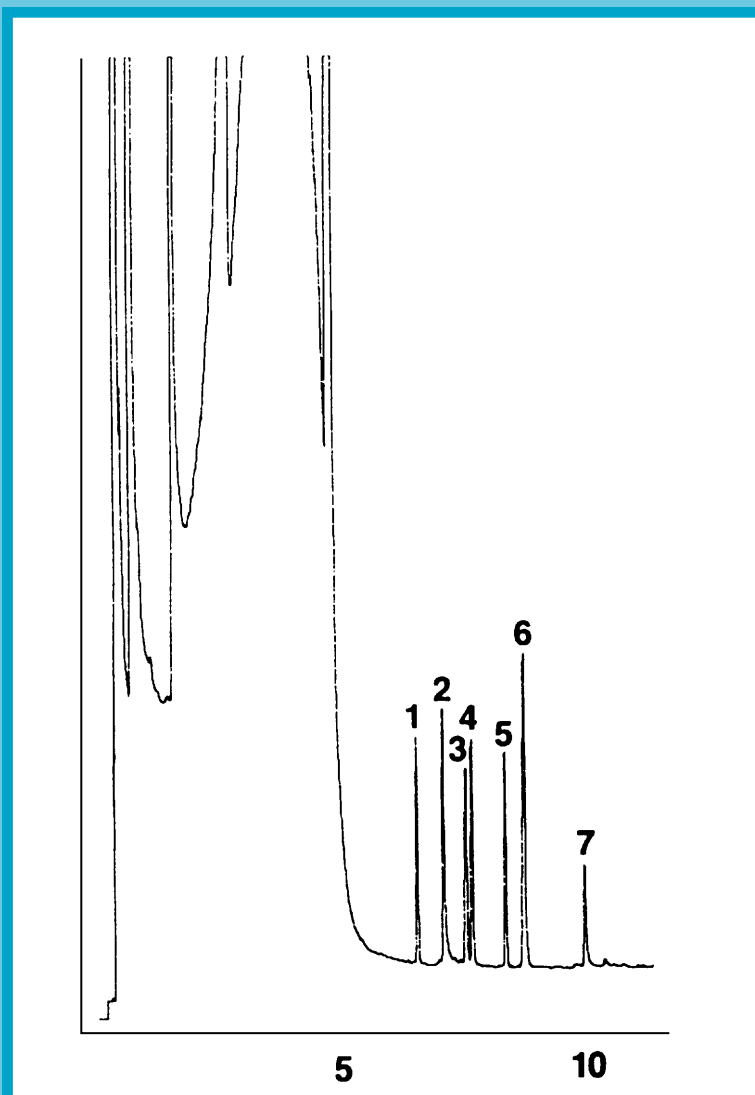

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com