



# GC CAPILLARY COLUMN

## APPLICATIONS: Ultra Alloy™ Stainless Steel

*Pesticides  
and  
Herbicides*

**COLUMN:** UAC-1, Dimethylpolysiloxane  
15M. x 0.25mm I.D. x 0.25µm film  
Cat. No.: UAC-1-15-0.25F


**Temperature:** 50° (1 min. hold) (20°/ min) -150° -  
(5°/min.) - 300°C

**Injector:** 250°C

**Detector:** 300°C, FID

**Carrier Gas:** 0.8 ml/min., Helium

- |                     |                    |                     |
|---------------------|--------------------|---------------------|
| 1. Chloneb          | 9. Terbucarb       | 17. Napropamide     |
| 2. Pencycuron       | 10. Fenitrothion   | 18. Butamifos       |
| 3. Benfluralin      | 11. Chlorpyrifos   | 19. Flutolanil      |
| 4. Simazine         | 12. Captan         | 20. Isoxathion      |
| 5. Chlorothalonil   | 13. Pendimethalin  | 21. Mepronil        |
| 6. Propyzamide      | 14. Methyl dymron  | 22. Pyridaphenthion |
| 7. Diazinon         | 15. Isofenphos     | 23. Iprodione       |
| 8. Tolclofos-methyl | 16. Isoprothiolane | 24. Bensulide       |

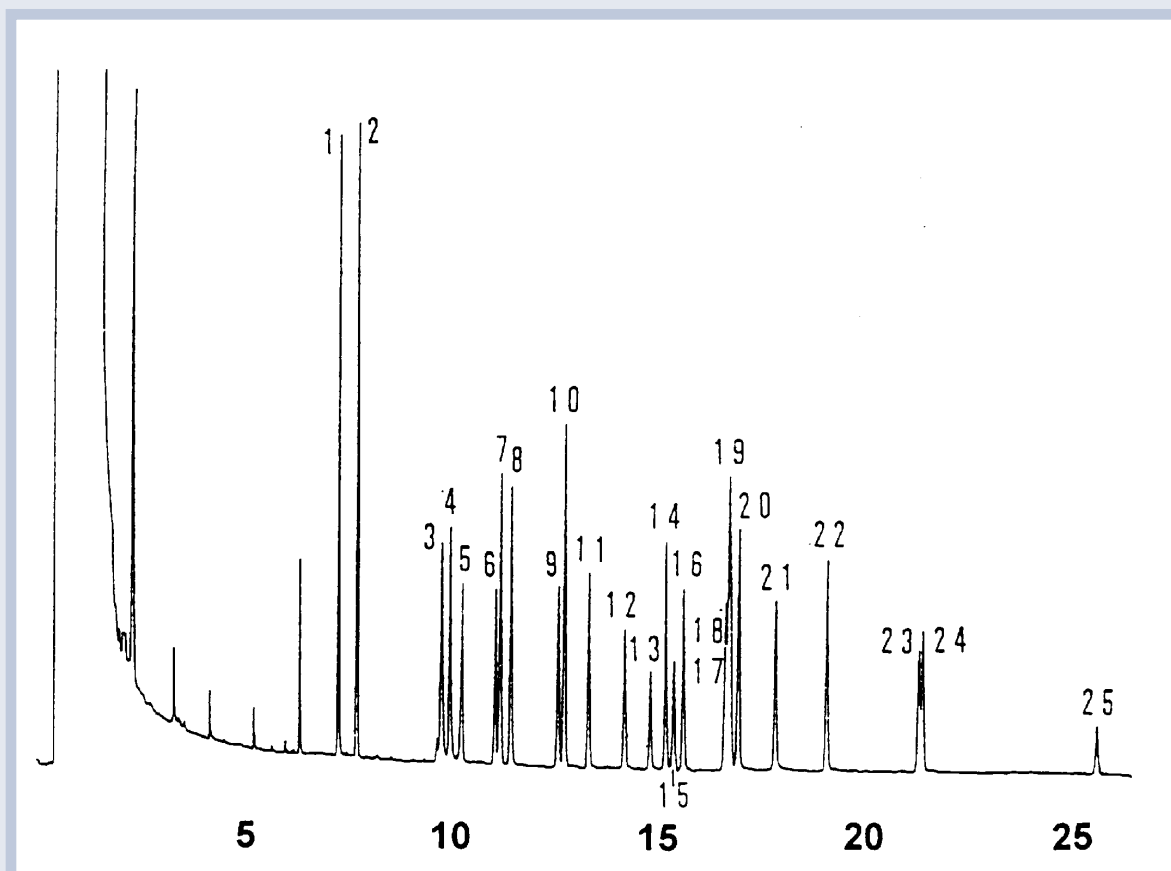

**QUADREX**



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com