



GC CAPILLARY COLUMN APPLICATIONS:

Petrochemical

#4 Fuel Oil

COLUMN: 007-5, (5% Phenyl) Methylpolysiloxane
25M. x 0.25mm I.D. x 0.25 μ m film
Cat. No.: 007-5-25-0.25F


Temperature: 80° (8°/min.) - 300°C

Injector: 250°C

Detector: 310°C, FID

Carrier Gas: 29 cm/sec., Helium

- | | |
|-------------|----------|
| 1. C-10 | 13. C-20 |
| 2. C-11 | 14. C-21 |
| 3. C-12 | 15. C-22 |
| 4. C-13 | 16. C-23 |
| 5. C-14 | 17. C-24 |
| 6. C-15 | 18. C-25 |
| 7. C-16 | 19. C-26 |
| 8. C-17 | 20. C-27 |
| 9. pristane | 21. C-28 |
| 10. C-18 | 22. C-29 |
| 11. phytane | 23. C-30 |
| 12. C-19 | |

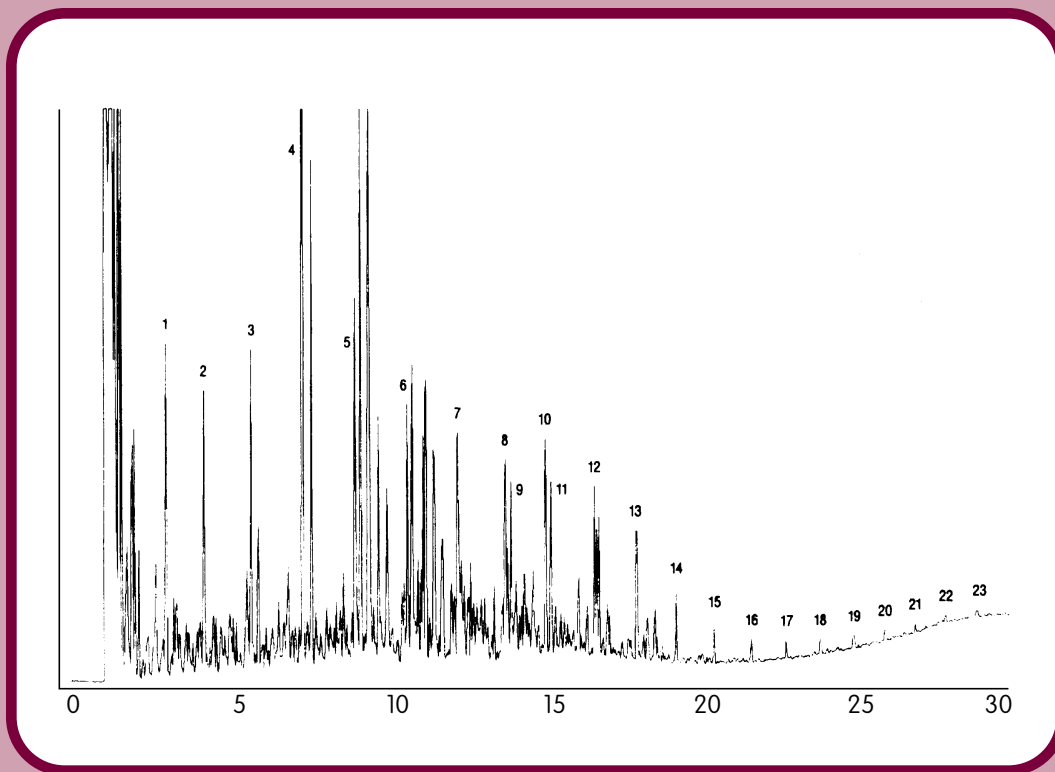

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com