




GC CAPILLARY COLUMN APPLICATIONS:

Petrochemical

Reference Gas Oil - ASTM 2887

COLUMN: 007-1, Dimethylpolysiloxane
10M. x 0.53mm I.D. X 1.0 μ m film
Cat. No.: 007-1-10V-1.0F

Temperature: 35° (25°/min.) - 320°C
Injector: 300°C
Detector: 350°C, FID
Carrier Gas: 20 ml/min., Helium

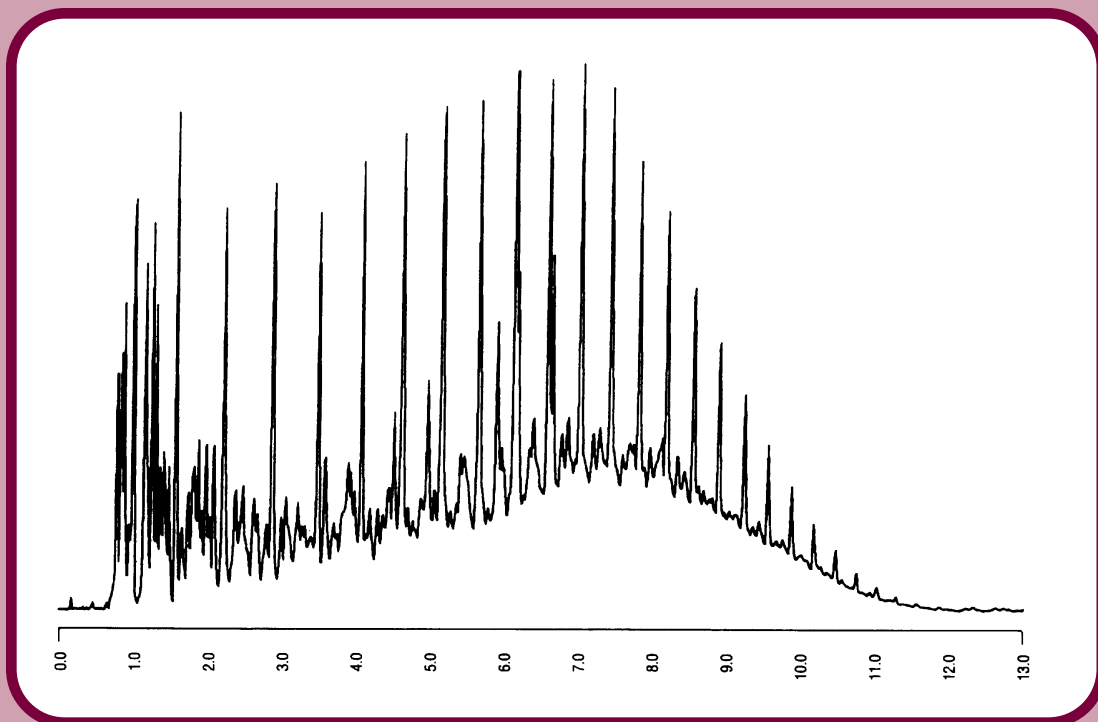

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com