




GC CAPILLARY COLUMN APPLICATIONS: Fames

Rapeseed Oil

COLUMN: 007-CW, Bonded Carbowax
25M. x 0.25mm I.D. x 0.25µm film
Cat. No.: 007-CW-25-0.25F

Temperature: 200° (2°/min.) - 240°C
Injector: 250°C
Detector: 250°C, FID
Carrier Gas: 28 cm/sec, Helium

1. C-14:0	4. C-18:1	7. C-20:0	10. C-22:1
2. C-16:0	5. C-18:2	8. C-20:1	11. C-24:0
3. C-18:0	6. C-18:3	9. C-22:0	

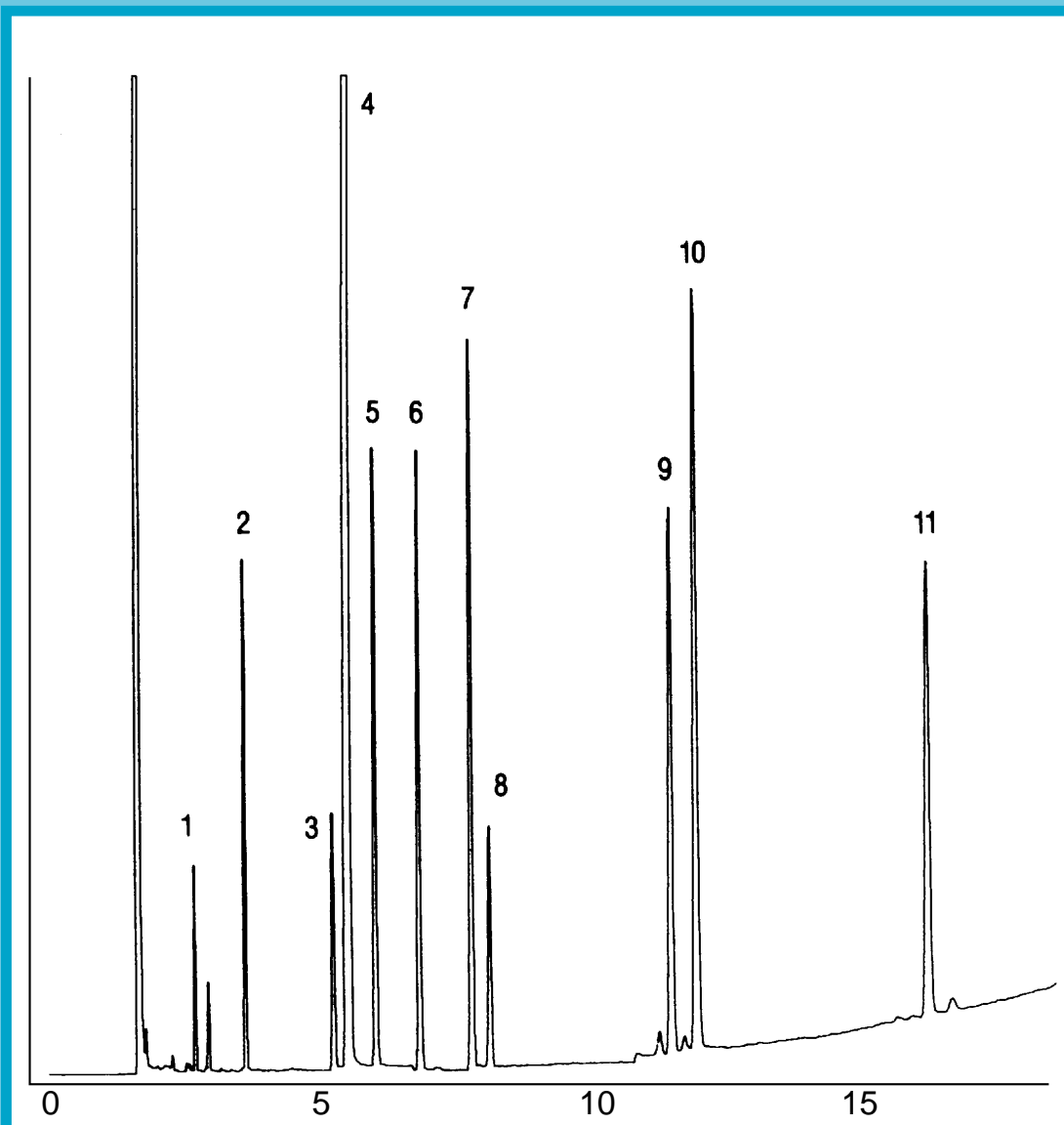

QUADREX



PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391
TOLLFREE: 800-275-7033 • E-MAIL: tech@quadrexcorp.com • WEB: www.quadrexcorp.com