




# GC CAPILLARY COLUMN APPLICATIONS: Fames

## Fame Mix

**COLUMN:** 007-CW, Bonded Carbowax  
25M. x 0.25mm I.D. x 0.25µm film  
Cat. No.: 007-CW-25-0.25F

**Temperature:** 80° (8°/min.) - 240°C  
**Injector:** 230°C  
**Detector:** 280°C, FID  
**Carrier Gas:** 28 cm/sec., Helium

1. C-8:0	6. C-14:0	11. C-18:1	16. C-20:2
2. C-9:0	7. C-14:1	12. C-18:2	17. C-20:3
3. C-10:0	8. C-16:0	13. C-18:3	
4. C-11:0	9. C-16:1	14. C-20:0	
5. C-12:0	10. C-18:0	15. C-20:1	

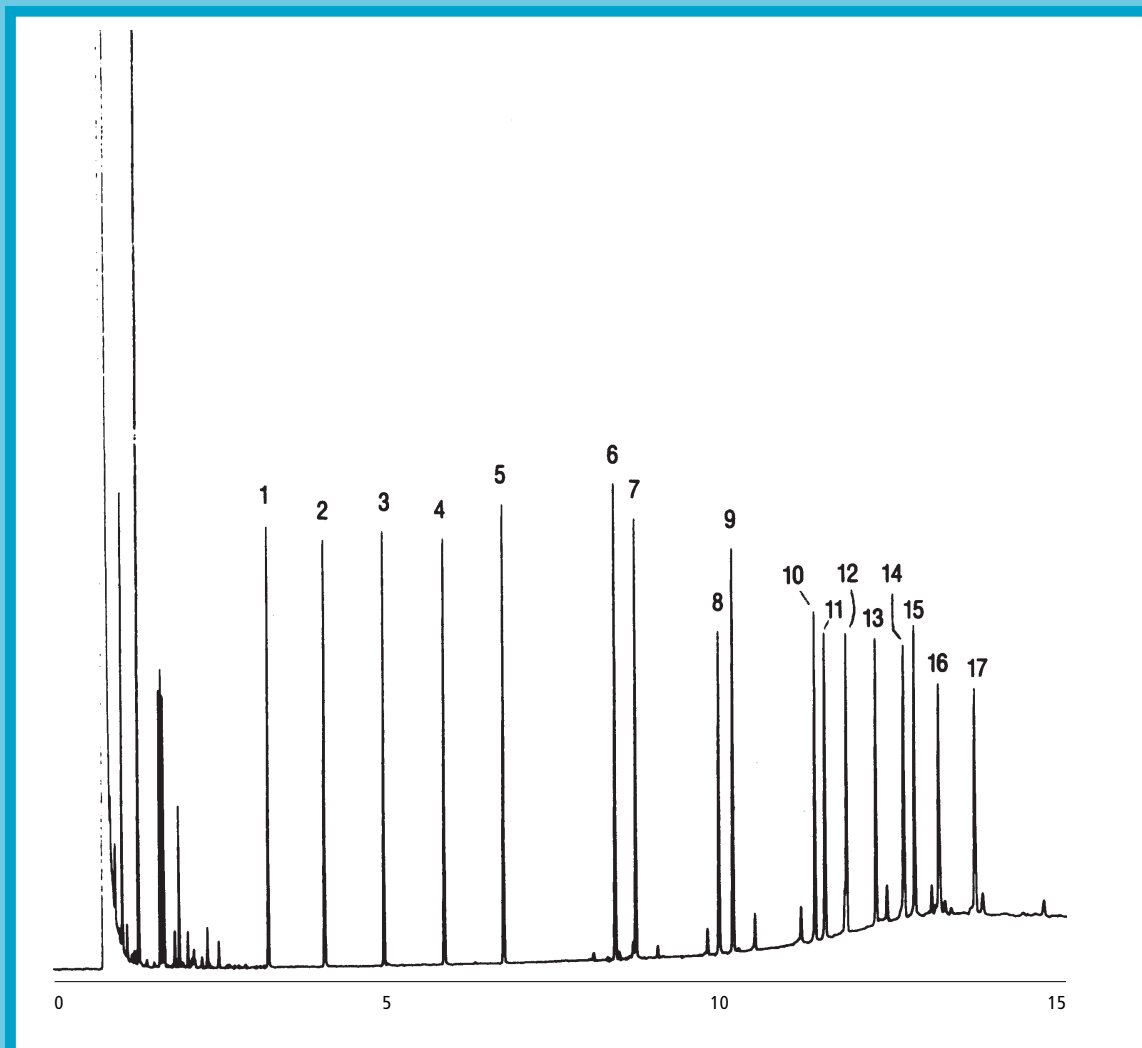

PO Box 3881, Woodbridge, CT 06525 • (203) 393-3112 • FAX (203) 393-0391  
TOLLFREE: 800-275-7033 • E-MAIL: [tech@quadrexcorp.com](mailto:tech@quadrexcorp.com) • WEB: [www.quadrexcorp.com](http://www.quadrexcorp.com)

**QUADREX**

